

Features & Benefits

High Quality, Detachable Universal Needle

- ▲ Provides optimum patient comfort
- ▲ Bevel-Up orientation for low angle injection
- ▲ Fits any conventional Luer Lock or Luer Slip syringe
- ▲ Allows for interchange of needles for aspiration and administration

ISO Color Coded Safety Cover

- ▲ Provides easy identification of Needle Gauge
- ▲ Two cover sizes to accommodate short and long needle lengths

Enhanced Safety Locking Mechanism

- ▲ Locking clips secure the needle inside the safety cover – two in the needle cover & two on the Secure-Lok™ mechanism
- ▲ Métier Medical's unique 'Secure-Lok™' mechanism provides enhanced locking at the base of the safety cover

Single-Handed Activation

- ▲ Requires minimal training
- ▲ Almost no change in hand position required to activate

Audible 'Click' Upon Activation

- ▲ Signals that cover is securely locked in place over the needle

Designed to safeguard the healthcare worker against the hazard of needlestick injury

Métier Medical™
Safe-T-Clip
Safety Hypodermic Needle

Providing
Simple
Solutions
for
Healthcare
Safety

 **Métier**
Medical Limited

Innovative Medical Safety Technologies

Lvl. 15, Corporate Centre One
2 Corporate Court, BUNDALL
QLD 4217 AUSTRALIA

Telephone: +61 7 5668 5436
Email: info@metiermedical.com

WWW.METIERMEDICAL.COM



Métier Medical™
Safe-T-Clip

Safety
Hypodermic
Needle




For use with any Conventional
Luer Lock or Luer Slip Syringe

Before Activation

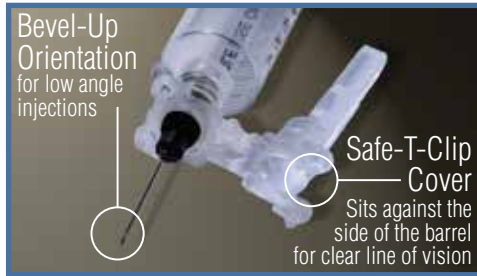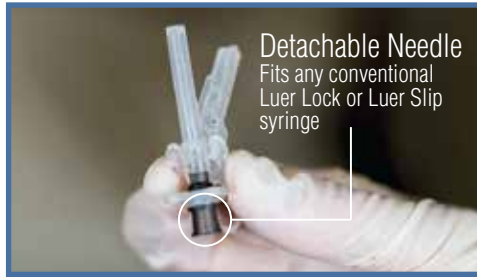

After Activation

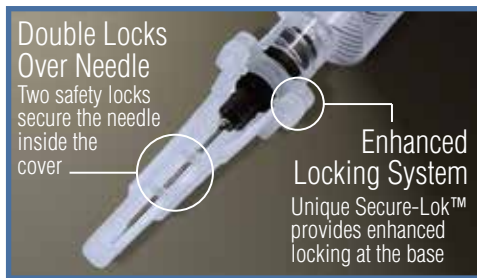

Métier Medical TM
Safe-T-Clip

Also available, Safe-T-Clip™ Syringe

The Safe-T-Clip™ needle is also supplied with a syringe option with safety needle mounted or bi-packed.

Safety Needle Product Options

PRODUCT CODE	SYRINGE SIZE	NEEDLE SIZE (Imperial)	NEEDLE SIZE (Metric)	PIECES PER BOX	PIECES PER CASE
Safe-T-Clip™ SAFETY HYPODERMIC NEEDLE (other Needle Size Variations Available)					
10212		30 G x 1/2"	0.30 mm x 12.7 mm	50	2500
10213		29 G x 1/2"	0.33 mm x 12.7 mm	50	2500
10378		27 G x 1"	0.4 mm x 25.4 mm	50	2500
10215		27 G x 1/2"	0.4 mm x 12.7 mm	50	2500
10216		27 G x 5/8"	0.4 mm x 15.87 mm	50	2500
10377		26 G x 1"	0.45 mm x 25.4 mm	50	2500
10217		26 G x 1/2"	0.45 mm x 12.7 mm	50	2500
10218		26 G x 5/8"	0.45 mm x 15.87 mm	50	2500
10219		25 G x 5/8"	0.5 mm x 15.87 mm	50	2500
10220		25 G x 3/4"	0.5 mm x 19.05 mm	50	2500
10221		25 G x 1"	0.5 mm x 25.4 mm	50	2500
10360		25 G x 1 1/2"	0.5 mm x 38.10 mm	50	2500
10376		24 G x 1 1/2"	0.55 mm x 38.10 mm	50	2500
10223		24 G x 1"	0.55 mm x 25.4 mm	50	2500
10375		23 G x 1 1/2"	0.60 mm x 38.10 mm	50	2500
10224		23 G x 1"	0.60 mm x 25.4 mm	50	2500
10361		22 G x 1 1/2"	0.70 mm x 38.10 mm	50	2500
10226		22 G x 1"	0.70 mm x 25.4 mm	50	2500
10362		21 G x 1 1/2"	0.80 mm x 38.10 mm	50	2500
10228		21 G x 1"	0.80 mm x 25.4 mm	50	2500
10374		20 G x 1 1/2"	0.90 mm x 38.10 mm	50	2500
10373		20 G x 1"	0.90 mm x 25.4 mm	50	2500
10372		19 G x 1 1/2"	1.10 mm x 38.10 mm	50	2500
10371		19 G x 1"	1.10 mm x 25.4 mm	50	2500
10370		18 G x 1 1/2"	1.20 mm x 38.10 mm	50	2500
10369		18 G x 1"	1.20 mm x 25.4 mm	50	2500
Safe-T-Clip™ NEEDLE WITH SYRINGE OPTION - LUER LOCK (with Mounted Safety Needle)					
10354	3 MI.	Syringe + Safety Needle		100	1200
10355	5 MI.	Syringe + Safety Needle		100	1200
10356	10 MI.	Syringe + Safety Needle		100	600
Safe-T-Clip™ NEEDLE WITH SYRINGE OPTION - LUER LOCK (with Bi-Packed Safety Needle)					
10357	3 MI.	Syringe + Safety Needle		100	1200
10358	5 MI.	Syringe + Safety Needle		100	1200
10359	10 MI.	Syringe + Safety Needle		100	600
Safe-T-Clip™ NEEDLE WITH SYRINGE OPTION - LUER SLIP (with Mounted Safety Needle)					
10363	3 MI.	Syringe + Safety Needle		100	1200
10364	5 MI.	Syringe + Safety Needle		100	1200
10365	10 MI.	Syringe + Safety Needle		100	600
Safe-T-Clip™ NEEDLE WITH SYRINGE OPTION - LUER SLIP (with Bi-Packed Safety Needle)					
10366	3 MI.	Syringe + Safety Needle		100	1200
10367	5 MI.	Syringe + Safety Needle		100	1200
10368	10 MI.	Syringe + Safety Needle		100	600

Instructions for Use

Step 01

Attach the Safe-T-Clip™ needle to any Luer Lock or Luer Slip syringe.



Step 02

Pull the safety cover back so that it sits parallel to the syringe barrel allowing for a clear view of the injection site. Remove the needle shield.



Step 03

Draw up medication in accordance with established protocol.



Step 04

Administer medication according to established protocol.



Step 05

Immediately following injection, the safety cover can be locked over the needle by depressing on a flat surface (an audible click indicates the safety mechanism has been activated). Discard in an approved sharps waste container.

